

The ***Ultimate***
Surgical Stapling Range

purple surgical

The *Ultimate* Reloadable Linear Cutter Stapler



The **Purple Surgical Ultimate** Reloadable Linear Cutter Stapler places two double staggered rows of titanium staples and simultaneously cuts and divides tissue between the two double rows.

The **Purple Surgical *Ultimate*** Reloadable Linear Cutter Stapler is available in 60mm and 80mm staple lengths with reloading units available for both standard tissue (3.8mm) and thick tissue (4.8mm). Each device may be reloaded up to 7 times for a maximum of 8 firings per patient. The ergonomic design, with rubber grip inlays and flexible dual lever for left or right handed staple firing, facilitates user comfort and control. Staple formation creates the perfect 'B' shape on each and every firing for improved clinical performance.

Features and benefits include:

- New blade with each cartridge guarantees smooth, precise transection with each firing, reducing the risk of cross infection
- Automatic safety lock-out mechanism incorporated within each cartridge
- Rear hinge allows for one handed operation
- Adjustable firing lever enables you to fire from either side of the device

Anatomy suitable for this device includes:

General Surgery - Stomach | Duodenum | Small Bowel | Liver | Colon
Gynaecological Surgery - Fallopian Tubes | Ovaries
Urologic Surgery - Bladder
Thoracic Surgery - Parenchyma

Ultimate Reloadable Linear Cutter Stapler

Product	Code	Colour	Tissue Type	Reload Included	Open Staple Height	Closed Staple Height	Pack Qty
Ultimate 60mm Reloadable Linear Cutter Stapler	PS9LC60B	Blue	Standard	PS9LCR60B	3.8mm	1.5mm	3 / Box
Ultimate 60mm Reloadable Linear Cutter Stapler	PS9LC60G	Green	Thick	PS9LCR60G	4.8mm	2.0mm	3 / Box
Ultimate 80mm Reloadable Linear Cutter Stapler	PS9LC80B	Blue	Standard	PS9LCR80B	3.8mm	1.5mm	3 / Box
Ultimate 80mm Reloadable Linear Cutter Stapler	PS9LC80G	Green	Thick	PS9LCR80G	4.8mm	2.0mm	3 / Box

Ultimate Linear Cutter Stapler Reloading Units

Product	Code	Colour	Tissue Type	Cut Length	Rows	Staples	Open Staple Height	Closed Staple Height	Pack Qty
Ultimate 60mm Linear Cutter Stapler Reloading Unit	PS9LCR60B	Blue	Standard	56mm	4	64	3.8mm	1.5mm	12 / Box
Ultimate 60mm Linear Cutter Stapler Reloading Unit	PS9LCR60G	Green	Thick	56mm	4	64	4.8mm	2.0mm	12 / Box
Ultimate 80mm Linear Cutter Stapler Reloading Unit	PS9LCR80B	Blue	Standard	76mm	4	84	3.8mm	1.5mm	12 / Box
Ultimate 80mm Linear Cutter Stapler Reloading Unit	PS9LCR80G	Green	Thick	76mm	4	84	4.8mm	2.0mm	12 / Box

The *Ultimate*

Reloadable Linear Stapler with Indicator



Cutting guide offers safety during transection

Low profile design for easy access

Ergonomic, single handed firing lever with rubber grips

Release button cancels firing process at stages 1, 2 or 3

Colour coded indication window:

- 1 Neutral position
- 2 Half closed position
- 3 Tissue compressed between cartridge & anvil
- 4 Staples fired

The **Purple Surgical Ultimate** Reloadable Linear Stapler with Indicator places two staggered rows of titanium staples, and is available in 30mm, 45mm, 60mm and 90mm staple line lengths with reload units available for both standard tissue (3.5mm) and thick tissue (4.8mm).

The **Purple Surgical *Ultimate*** Reloadable Linear Stapler with Indicator may be reloaded up to 7 times for a maximum of 8 firings per patient. The ergonomic design, with rubber grip and single handed firing lever, facilitates user comfort and control. Staple formation creates the perfect 'B' shape on each and every firing for improved clinical performance.

Features and benefits include:

- Cutting guide for enhanced safety during transection
- Colour coded indication window for easy cancelation at stages 1 to 3
- Retaining pin can either be manually or mechanically engaged
- Available in 4 size options (30mm, 45mm, 60mm and 90mm)

Anatomy suitable for this device includes:

General Surgery - Oesophagus | Stomach | Duodenum | Small Bowel | Liver | Colon | Pylorus
Gynaecological Surgery - Fallopian Tubes | Ovaries
Thoracic Surgery - Parenchyma | Bronchus

Ultimate Reloadable Linear Stapler with Indicator

Product	Code	Colour	Tissue Type	Reload Included	Open Staple Height	Closed Staple Height	Pack Qty
<i>Ultimate</i> 30mm Reloadable Linear Stapler with Indicator	PS9LSI30B	● Blue	Standard	PS9LSR30B	3.5mm	1.5mm	3 / Box
<i>Ultimate</i> 30mm Reloadable Linear Stapler with Indicator	PS9LSI30G	● Green	Thick	PS9LSR30G	4.8mm	2.0mm	3 / Box
<i>Ultimate</i> 45mm Reloadable Linear Stapler with Indicator	PS9LSI45B	● Blue	Standard	PS9LSR45B	3.5mm	1.5mm	3 / Box
<i>Ultimate</i> 45mm Reloadable Linear Stapler with Indicator	PS9LSI45G	● Green	Thick	PS9LSR45G	4.8mm	2.0mm	3 / Box
<i>Ultimate</i> 60mm Reloadable Linear Stapler with Indicator	PS9LSI60B	● Blue	Standard	PS9LSR60B	3.5mm	1.5mm	3 / Box
<i>Ultimate</i> 60mm Reloadable Linear Stapler with Indicator	PS9LSI60G	● Green	Thick	PS9LSR60G	4.8mm	2.0mm	3 / Box
<i>Ultimate</i> 90mm Reloadable Linear Stapler with Indicator	PS9LSI90B	● Blue	Standard	PS9LSR90B	3.5mm	1.5mm	3 / Box
<i>Ultimate</i> 90mm Reloadable Linear Stapler with Indicator	PS9LSI90G	● Green	Thick	PS9LSR90G	4.8mm	2.0mm	3 / Box

Ultimate Linear Stapler Reloading Units

Product	Code	Colour	Tissue Type	Rows	Staples	Open Staple Height	Closed Staple Height	Pack Qty
<i>Ultimate</i> 30mm Linear Stapler Reloading Unit	PS9LSR30B	● Blue	Standard	2	11	3.5mm	1.5mm	12 / Box
<i>Ultimate</i> 30mm Linear Stapler Reloading Unit	PS9LSR30G	● Green	Thick	2	11	4.8mm	2.0mm	12 / Box
<i>Ultimate</i> 45mm Linear Stapler Reloading Unit	PS9LSR45B	● Blue	Standard	2	15	3.5mm	1.5mm	12 / Box
<i>Ultimate</i> 45mm Linear Stapler Reloading Unit	PS9LSR45G	● Green	Thick	2	15	4.8mm	2.0mm	12 / Box
<i>Ultimate</i> 60mm Linear Stapler Reloading Unit	PS9LSR60B	● Blue	Standard	2	21	3.5mm	1.5mm	12 / Box
<i>Ultimate</i> 60mm Linear Stapler Reloading Unit	PS9LSR60G	● Green	Thick	2	21	4.8mm	2.0mm	12 / Box
<i>Ultimate</i> 90mm Linear Stapler Reloading Unit	PS9LSR90B	● Blue	Standard	2	33	3.5mm	1.5mm	12 / Box
<i>Ultimate</i> 90mm Linear Stapler Reloading Unit	PS9LSR90G	● Green	Thick	2	33	4.8mm	2.0mm	12 / Box

The *Ultimate* Circular Stapler with Dual Safety



The **Purple Surgical *Ultimate*** Circular Stapler with Dual Safety places a circular staggered double row of titanium staples in the tissue and resects excess tissue to create a circular anastomosis.

The **Purple Surgical *Ultimate*** Circular Stapler with Dual Safety is available in 5 size options, indicated by colour coded anvil heads. 3 staple sizes are available to provide optimal tissue thickness matching. The ergonomic design and textured firing handle facilitates user comfort and control. Staple formation creates the perfect 'B' shape staple when the red line is positioned centrally within the green firing indicator window.

Features and benefits include:

- Dual safety lock-out mechanism prevents misfiring
- Atraumatic Dilating Tip Trocar for reduced trauma
- Fin design Wing Nut for easy counting of full or half turns
- Available in 5 size options with 3 open staple heights and variable closed staple heights for greater choice and flexibility

Anatomy suitable for this device includes:

General Surgery • Oesophagus | Stomach | Duodenum | Colon | Rectum



Ultimate Circular Stapler with Dual Safety

Product	Code	Colour	Lumen	Shaft Length	Rows	Staples	Open Staple Height	Closed Staple Height	Pack Qty
<i>Ultimate</i> 24mm Circular Stapler with Dual Safety	PS9CS24S	● Yellow	15mm	18.5cm	2	18	4.5mm	1.8mm	3 / Box
<i>Ultimate</i> 26mm Circular Stapler with Dual Safety	PS9CS26S	○ White	17mm	18.5cm	2	20	4.8mm	2.0mm	3 / Box
<i>Ultimate</i> 29mm Circular Stapler with Dual Safety	PS9CS29S	● Blue	20mm	18.5cm	2	24	4.8mm	2.0mm	3 / Box
<i>Ultimate</i> 32mm Circular Stapler with Dual Safety	PS9CS32S	● Green	24mm	18.5cm	2	30	5.0mm	2.2mm	3 / Box
<i>Ultimate</i> 34mm Circular Stapler with Dual Safety	PS9CS34S	● Purple	24.5mm	18.5cm	2	32	5.0mm	2.2mm	3 / Box

The *Ultimate* Circular Stapler with Dual Safety & Tilt Top Anvil



Low profile Tilt
Top anvil for easy
insertion and
removal

Green firing
indicator window
offering variable
closed staple
heights

Angled quick release
wing nut for faster
opening and closing

Dual safety
lock-out mechanism
prevents misfiring

Vented anvil for
vacuum release

The **Purple Surgical *Ultimate*** Circular Stapler with Dual Safety & Tilt Top Anvil places a circular staggered, double row of titanium staples in the tissue and resects excess tissue to create a circular anastomosis.

The **Purple Surgical Ultimate** Circular Stapler with Dual Safety & Tilt Top anvil is available in 4 size options indicated by colour coded anvil washers. 3 staple sizes are available to provide optimal tissue thickness matching. The ergonomic design and textured firing handle facilitates user comfort and control. Staple formation creates the perfect 'B' shape staple when the red line is positioned centrally within the green firing indicator window.

Features and benefits include:

- Dual safety lock-out mechanism prevents misfiring
- Vented, low profile tilt top anvil releases vacuum and minimises the potential of tissue trauma on removal
- Angled wing nut for easy counting of full and half turns
- Available in 4 size options with 3 open staple heights and variable closed staple heights for greater choice and flexibility

Anatomy suitable for this device include:


General Surgery | Oesophagus | Stomach | Duodenum | Colon | Rectum



Ultimate Circular Stapler with Dual Safety & Tilt Top anvil

Product	Code	Colour	Lumen	Shaft Length	Rows	Staples	Open Staple Height	Closed Staple Height	Pack Qty
Ultimate 24mm Circular Stapler with Dual Safety & Tilt Top Anvil	PS9CSA24S	● Yellow	15mm	18.5cm	2	18	4.5mm	1.8mm	3 / Box
Ultimate 26mm Circular Stapler with Dual Safety & Tilt Top Anvil	PS9CSA26S	○ White	17mm	18.5cm	2	20	4.8mm	2.0mm	3 / Box
Ultimate 29mm Circular Stapler with Dual Safety & Tilt Top Anvil	PS9CSA29S	● Blue	20mm	18.5cm	2	24	4.8mm	2.0mm	3 / Box
Ultimate 32mm Circular Stapler with Dual Safety & Tilt Top Anvil	PS9CSA32S	● Green	24mm	18.5cm	2	30	5.0mm	2.2mm	3 / Box

The *Ultimate* Haemorrhoidal Circular Stapler with Dual Safety & Transparent Cartridge

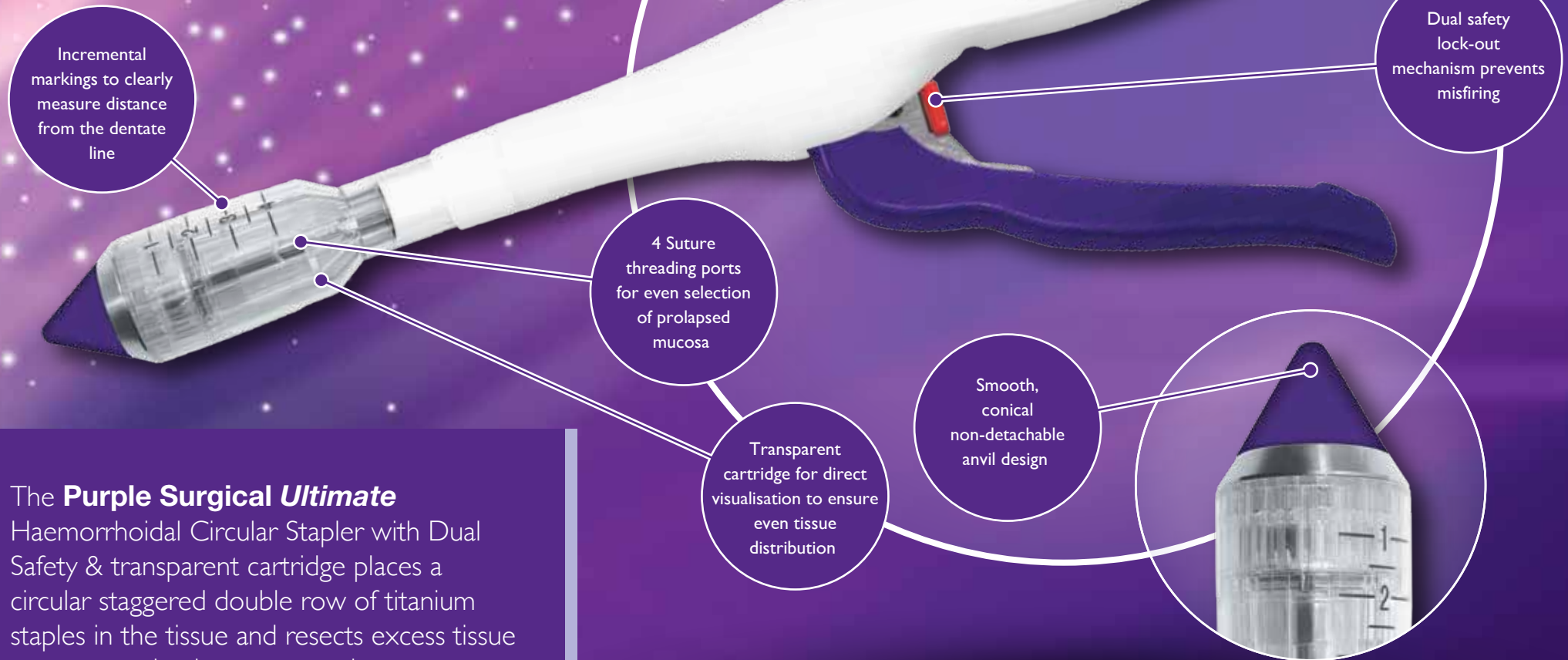

Green firing
indicator
window

Dual safety
lock-out
mechanism prevents
misfiring

Incremental
markings to clearly
measure distance
from the dentate
line

4 Suture
threading ports
for even selection
of prolapsed
mucosa

Transparent
cartridge for direct
visualisation to ensure
even tissue
distribution

Smooth,
conical
non-detachable
anvil design



The **Purple Surgical Ultimate**
Haemorrhoidal Circular Stapler with Dual
Safety & transparent cartridge places a
circular staggered double row of titanium
staples in the tissue and resects excess tissue
to create a circular anastomosis.

The **Purple Surgical *Ultimate*** Haemorrhoidal Circular Stapler with Dual Safety & transparent cartridge is designed for use in the surgical treatment of haemorrhoidal disease. Available in two diameters; 32/33mm and 33/34mm and supplied with 5 accessories including both a clear circular anal dilator and clear butterfly anal dilator.

Features and benefits include:

- Dual safety lock out mechanism prevents misfiring
- Fixed anvil to eliminate risk of detachment from the cartridge
- Transparent cartridge for direct visualisation of tissue within the cartridge
- Double staggered row of staples, for secure staple line formation

The ***Ultimate*** Haemorrhoidal Circular Stapler is supplied with the following 5 accessories:



Clear Circular Anal Dilator



Clear Butterfly Anal Dilator



Purse-String Suture Anoscope



Anal Dilator Obturator



Suture Threader

Ultimate Haemorrhoidal Circular Stapler with Dual Safety & Transparent Cartridge

Product	Code	Rows	Staples	Knife Diameter	Open Staple Height	Closed Staple Height	Pack Qty
<i>Ultimate</i> 32/33mm Haemorrhoidal Circular Stapler	PS9HCS32/33	2	32	22mm	3.8mm	1.5mm	3 / Box
<i>Ultimate</i> 33/34mm Haemorrhoidal Circular Stapler	PS9HCS33/34	2	34	24mm	3.8mm	1.5mm	3 / Box



Purple Surgical International Limited | Purple Surgical UK Limited
2 Chestnut House, Farm Close, Shenley, Herts WD7 9AD, England
Tel: +44 (0)1923 839333 Fax: +44 (0)1923 839444 Email: mail@purplesurgical.com

www.purplesurgical.com